

# accutest® TSH

## A new test for Hypothyroidism that can be done in your office...





**Almost 27 million Americans are hypothyroid... less than 50% have been diagnosed!** Symptoms of Hypothyroidism are common to many clinical conditions and include:

- Weight gain, constipation
- Depression, lethargy, fatigue
- Sluggishness, sleepiness
- Menstrual irregularities
- Muscle cramps and aches
- Cold intolerance




**If you suspect hypothyroidism, test with CLIA-waived Accutest TSH.**

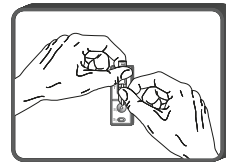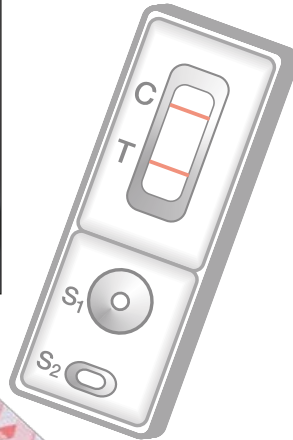

Place one drop of blood from a finger-stick or venous draw in S1



Wait 90 seconds

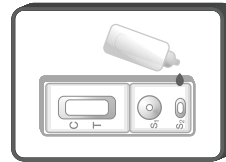

Place 4 drops of buffer in S2



Read results in 10 minutes




**American Thyroid Association recommends TSH as initial screen for hypothyroidism!**

- **CLIA - Waived** • **Cost Effective**
- **Easy to Use:** One drop of blood from finger-stick or venous blood draw
- **Rapid Results:** Know in 10 minutes or less if TSH is elevated
- **Accurate:** 98.5% Sensitivity and 96.9% Specificity
- **Room Temperature Storage**
- **Catalog No. MD801** • **CPT Code: 84443QW**
- **Reimbursement: \$23.47 in most States.**  
**Some private insurances higher.**



DISTRIBUTED BY:

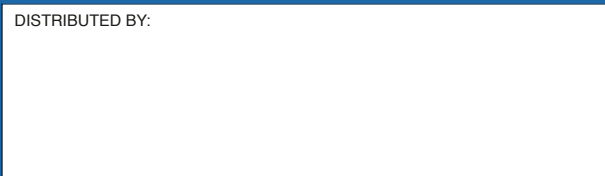

**JANT PHARMACAL CORPORATION**  
16255 Ventura Blvd., #505, Encino, CA 91436  
800.676.5565 818.986.8530 Fax 818.986.0235  
www.accutest.net info@accutest.net

Accutest is a registered trademark of Jant Pharmacal Corporation